






The innovative, patented device designed for greater volume collection that enables higher-sensitivity testing and eliminates the need for cold chain.

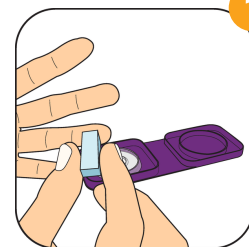


Stick finger with lancet.

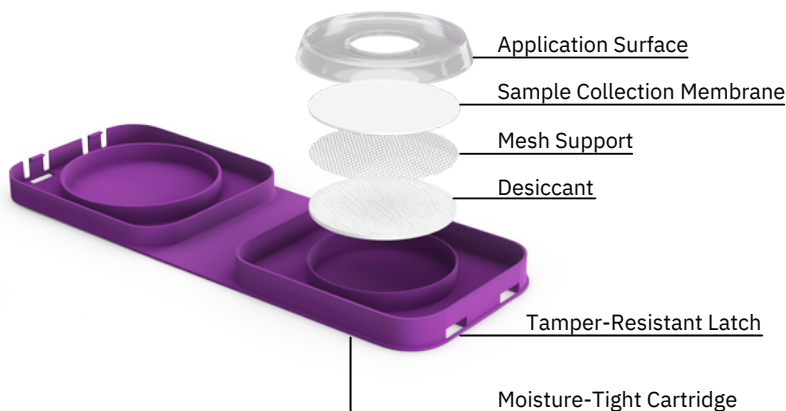

Add 5 drops of blood.

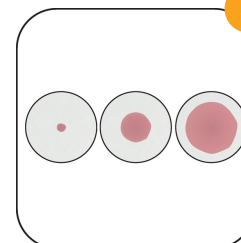


Wait 1 minute.

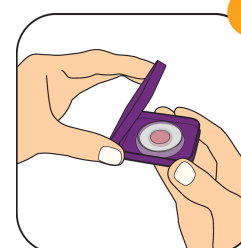

Close lid and ship.

Unique Collection Membrane


HemaSpot™ HD is designed for high volume sample collection, up to 160 µL. The larger sample collection membrane allows for multiple punches, providing the opportunity to run a variety of assays on a single sample.

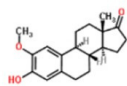


Small Molecules
Hormones, lipids, metabolites
LC-MS/MS, HPLC
Drug levels



Nucleic acids
DNA, RNA, mRNA, miRNA
RT-PCR, NGS
Viral load

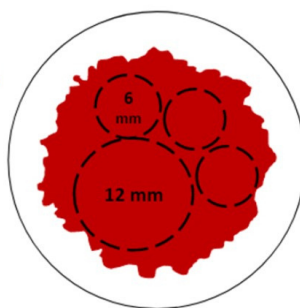

Proteins
Cytokines, enzymes, antibodies
ELISA, XMAP, MS
Serology

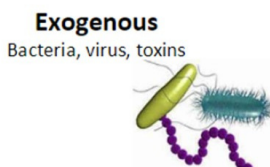

Exogenous
Bacteria, virus, toxins

For Research Use Only. Not For Use in Diagnostic Procedures.



LARGER MEMBRANE

Holds up to 160µL of bodily fluids; a greater volume than alternatives.



HIGH SENSITIVITY

Ideal for tests that require a greater volume of biospecimens.



MULTIPLE PUNCHES

Multiple tests can be performed with a single sample.



FLEXIBILITY

Suited for a wide variety of liquid biospecimens.