





The easy to use, patented device for collection and transportation of microsamples that can be safely stored at ambient temperatures.

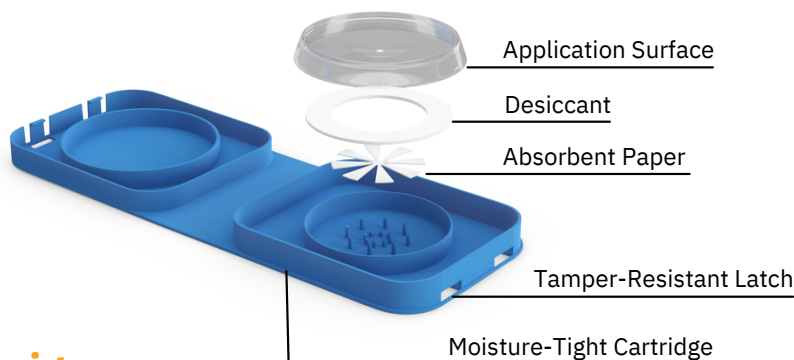

## Unique Collection Design


HemaSpot HF contains the innovative HemaForm™ absorbent paper. The unique fan-shaped design promotes even distribution across each of the eight blades allowing for reproducible analysis from the sample – no punching required.

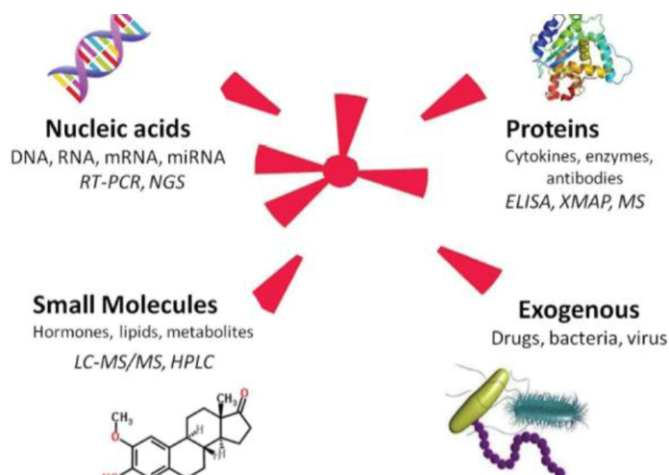

*Perform a variety of downstream assays at different points in time on a single sample.*

*Archive unused blades for future investigation.*


**For Research Use Only. Not For Use in Diagnostic Procedures.**

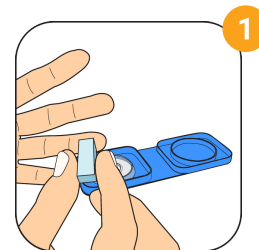


Stick finger with lancet.



Add 2 drops of blood.

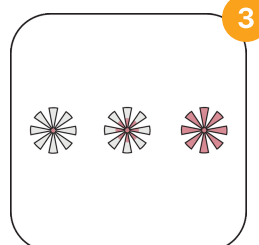

Wait 1 minute.

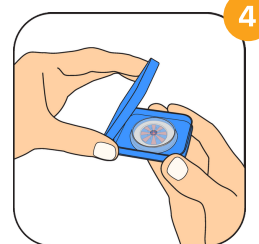

Close lid and ship.



### EVEN DISTRIBUTION

Fluid disperses uniformly across blades, reducing hematocrit effect.



### REPLICATE SAMPLES

Eight identical blades allow for repetitive and reproducible testing.



### PUNCHING NOT REQUIRED

Reduces possibility for contamination during process.



### NO SAMPLE WASTE

Pre-cut fan form enables 100% sample utilization.